






STANDING HEIGHT

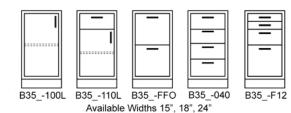

Available Widths 30", 36", 42", 48"

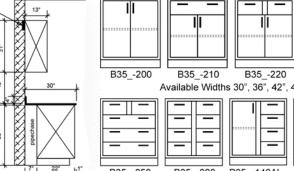

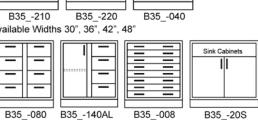

Also 52", 58"

CORNER

Cabinets are available in

modular steel, high quality

wood, stainless steel and

plastic laminate.


CHEMICAL

ACID STORAGE ASB35 -200 FSB35_-200




WALL



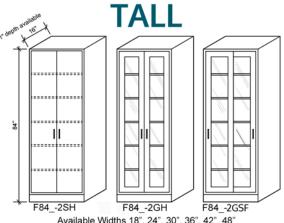

Available Widths 18", 24", 30", 36", 42", 48" (Note: 18" & 24" have one door)

CABINET DETAILS

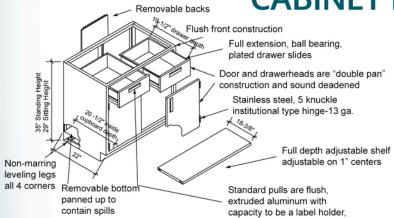

All surfaces shall be painted or zinc plated. The paint is chemical resistant baked-on epoxy powder. Specific test data is available upon request.

Gauges of cabinet components are as follows:

Aprons	18 ga.	Legs	18 ga.
Back Panels	20 ga.	Shelves	18 ga.
Bottom Panels	18 ga.	Side Panels	18 ga.
Drawer/Door Outer Panel	20 ga.	Table Frames	18 ga.
Door Inner Panel	20 ga.	K & V Shelf Clip	р
Drawer Bodies	20 ga.		


For complete specifications, request our full line of product specifications. Subject to change.




TABLES




Tubular Steel Tables - 16 ga. fully welded tables available in fixed or adjustable height with leveling legs or casters

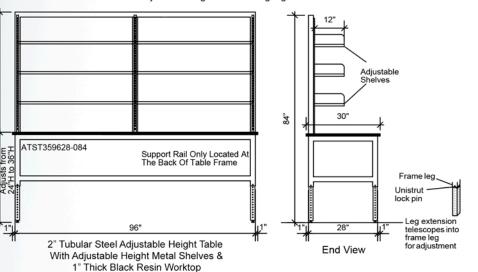

Tables are available in welded tubular steel, modular steel, stainless steel, and wood.

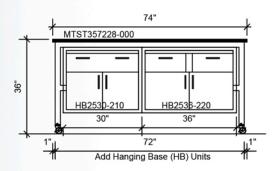


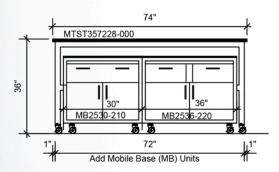

LEGS/APRONS



Combine aprons and legs to make the table of your size UT354822-002 = 35"h x 48" W x 22"d with 2 x drawers

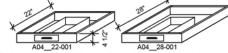

Available Widths 24", 30", 36", 42", 48", 54", 60", 66", 72"


COUNTERTOP SELECTION AND COMPARISON

3322						-	-	
	Chemical Resistance	Abrasion Resistance	Heat Resistance	Impact Resistance	Moisture Resistance	Bacterial Resistance	Areas of Use	Cost
Epoxy Resin	Excellent	Good	Excellent	Good	Excellent	Excellent	All	Medium
Phenolic Resin	Excellent	Good	Good	Good	Excellent	Excellent	All	Medium
Stainless Steel	Good	Good	Good	Excellent	Excellent	Excellent	Select	High
Edge Grain Maple	Fair	Fair	Fair	Good	Fair	Poor	Select	Medium
High Pressure Laminate	Fair	Good	Fair	Good	Fair	Fair	Select	Low






FUME HOODS

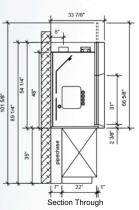


+	12" dia rour outlet coi exhaust blower ; & ductwork by others	
	B-606 - 72" Airfoil Bypass Fume Hood	8
	F-200 Light SW F-100 Gas F-100 Air X F-100 Vac X F-100 N2	F-278 Alarm V Blogg Oye St I F-100 Alf F-100 Vac × F-100 N X F-100 W
	F-225 elec	F 225 elec FLAMMABLE KEEP FIRE AWAY
	ASB3536-200	FSB3536-200

Mechanica	al Specificatiions	Sas	h Full Open @ 100 FPM	9	
Model	Hood Width	Exhaust Collar	CFM	S.P.	Sash Opening
B-604	48"	9" Round	867	.39"	31" x 38"
B-605	60"	10" Round	1,133	.38"	31" x 50"
B-606	72"	12" Round	1,400	.37"	31" x 62"
B-607	84"	2 x 9" Round	1,667	.37"	31' x 74"
B-608	96"	2 x 9" Round	1,933	.45"	31" x 86"

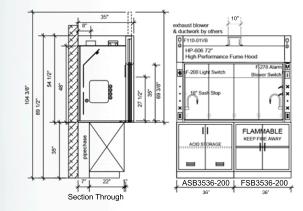

Hood Width	Exhaust Collar	CFM	S.P.	Sash Opening
48"	8" Round	380	.13"	27-1/2" x 38"
60"	8" Round	500	.19"	27-1/2" x 50"
72"	10" Round	620	.14"	27-1/2" x 62"
84"	2 x 8" Round	740	.12"	27-1/2' x 74"
96"	2 x 8" Round	860	.15"	27-1/2" x 86"
	48" 60" 72" 84"	48" 8" Round 60" 8" Round 72" 10" Round 84" 2 x 8" Round	48" 8" Round 380 60" 8" Round 500 72" 10" Round 620 84" 2 x 8" Round 740	48" 8" Round 380 .13" 60" 8" Round 500 .19" 72" 10" Round 620 .14" 84" 2 x 8" Round 740 .12"

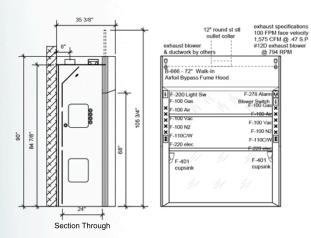

Mechanical Specifications		One	pen @		
Model	Hood Width	Exhaust Collar	CFM	S.P.	Sash Opening
B-9660	72"	12" Round	1,575	.47"	31" x 62"
B-9680	96"	2 x 9" Round	2 175	57"	31" x 86"

AIRFOIL BYPASS HOODS

Airfoil Bypass Hoods have the following Standard Items:

- 1. White Resin-Chem lining.
- 2. 16 ga. internal structural frame.
- 3. Non-corrosive nylon fasteners no metal inside the hood.
- 4. Full view 4-sided PVC vertical sash with full width pull and 7/32" laminated safety glass panel.
- 5. Optional manual sash lock at 18" operational height. 6. Lower airfoil with black epoxy finish.
- 7. Vapor proof T-8 fluorescent light fixture, with bulbs.
- 8. Round stainless steel duct collar.
- 9. Pre-punched fixture holes in post and inner side liner.
- 10. Gasketed access panel.
- 11. Black PVC sash guides.
- 12. Pre-wired electric fixtures to J-box on top of hood.
- 13. Pre-plumb service fixtures up or down verify.
- 14. UL-1805 Classified.

HIGH PERFORMANCE HOODS

Consumes 50% less energy than conventional hoods

HP Hoods have the following Standard Items:

- 1. White Resin-Chem lining.
- 2. 16 ga. internal structural frame.
- 3. Non-corrosive nylon fasteners no metal inside the hood.
- 4. Full view combination horizontal / vertical sash with top hung 7/32" laminated safety glass panels.
- 5. Keyed sash lock at 18" operational height.
- 6. Combination airfoil / flush sill containment trough with white plas-cote finish.
- 7. Vapor proof T-8 fluorescent light fixture, with bulbs.
- 8. Round stainless steel duct collar.
- 9. Pre-punched fixture holes in post and inner side liner.
- 10. Gasketed access panel.
- 11. Black PVC sash guides.
- 12. Pre-wired electric fixtures to J-box on top of hood.
- 13. Pre-plumb service fixtures up or down verify.
- 14. UL-1805 Classified.


WALK-IN AIRFOIL BYPASS HOODS

Walk-In Airfoil Bypass Hoods have the following Standard Items:


- 1. White Resin-Chem lining.
- 2. 16 ga. internal structural frame.
- 3. Non-corrosive nylon fasteners no metal inside the hood.
- 4. Full view 4-sided PVC vertical sash with full width pull and 7/32" laminated safety glass panel.
- 5. Optional manual sash lock at 18" operational height.
- 6. Lower airfoil with black epoxy finish.
- 7. Vapor proof T-8 fluorescent light fixture, with bulbs.
- 8. Round stainless steel duct collar.
- 9. Pre-punched fixture holes in post and inner side liner.
- 10. Gasketed access panel.
- 11. Black PVC sash guides.
- 12. Pre-wired electric fixtures to J-box on top of hood.
- 13. Pre-plumb service fixtures up or down verify.
- 14. UL-1805 Classified.

ACCESSORIES


EPOXY RESIN SINKS

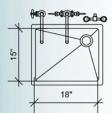

PL-30 Lipped 18" x 15" x 11" PSO-3 outlet

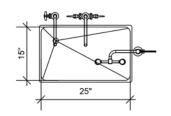


PL-55 Lipped 25" x 15" x 10" PSO-3 outlet

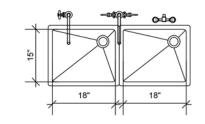


2 x PL-30 Lipped 18" x 15" x 11" / each 2 x PSO-3 outlets A large variety of sinks, fixtures, and safety items are available.


STAINLESS STEEL SINKS

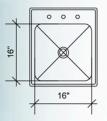

SL-2219-A-GR Ledge Type 16" x 16" x 7-1/2" J-35 outlet

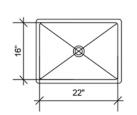

SX-1925-A No-Ledge 16" x 22" x 10-1/2" J-35 outlet

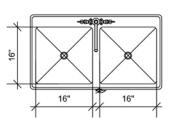

DL-2137-A-GR Double Bowl Ledge Type 16" x 16" x 7-1/2" / each 2 x J-35 outlets


SAFETY EQUIPMENT




EW-1022 Deck Mtd. Eyewash/Drench Hose



EW-806 Deck Mtd. Swivel Eyewash

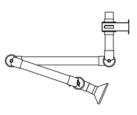

SH-75 Point of Use Snorkel Hood

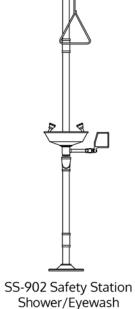

SINK FIXTURES



L-611 Cold Water



L-7611 DI Pure Water



L-412-VB-BH-08 H/C Water



L-512-WS Panel Mtd. H/C Water



L-2212-VB H/C Water



ELECTRIC/GAS FIXTURES



E-533 Double Gang Single Sided



L-4200-132S Double 180°



L-4200-131 Single





Available in Various Sizes & Materials (stainless steel, resin, acrylic)



LABORATORY FURNITURE LLC

Phone: 262-204-7660 Fax: 262-661-7060 www.psalaboratoryfurniture.com



DESIGNING LABORATORIES WITH YOU IN MIND AT HOME AND AROUND THE WORLD!

The tables and cabinetry met / exceeded our expectations in regards to quality and overall "fit" in the lab, and the setup process went very smoothly overall. We will definitely keep your team in mind for future projects.

- Aldevron

Thank you again for the wonderful work that was done. You guys have a great team and I will be happy to recommend you to anyone looking for a great furnishing company!

- Ayass Lung Clinic

We are very pleased with the custom new cabinets and benches in our MSI production area. Everything was engineered to fit perfectly. PSA did an excellent job and it was a pleasure to work with PSA on the installation. The quality of construction is exactly what I was expecting and hoping for. I am really impressed with the black phenolic resin tops, this material is the real deal when it comes to classic lab benchtops.

- Saint-Gobain

Agricultural

- -Pan American Seeds, Elburn, IL
- -ABS Global, De Forest, WI
- -Arkansas Mill Supply, White Hall, AR
- -Aldevron, Fargo, ND
- -Phibro Animal Health, Quincy, IL
- -Monsanto, Ankeny, IA

Research & Bio-Tech

- -Novozymes BioAg Inc., Milwaukee, WI
- -Synata Bio, Warrenville, IL
- -Virent, Madison, WI
- -Kurion Research Center, Richland, WA
- -Thermo Fisher Scientific

Pharmaceutical

- -Standard Process Laboratories, Palmyra, WI
- -Carma Labs, Franklin, WI
- -Agnitio (Somnium Life Sciences), Appleton, WI
- -Deerland Enzymes, Inc., Kennesaw, GA
- -STI Pharma LLC, Newtown, PA
- -Baxter HealthCare, Round Lake, IL
- -Catalent Pharmaceuticals, Madison, WI

Food & Beverage

- -California Dairies, Visalia, CA
- -Tate & Lyle, US Nationwide
- -Archer Daniels Midland Company, US Nationwide
- -Robert Mondovi Winery, Napa, CA
- -PepsiCo, City of Industry, CA
- -Tropicana Foods, FL
- -Ajinomoto Heartland, Eddyville, IA
- -Sweeteners Supply, Brookfield, IL
- -Golden Peanut, Dawson, GA
- -PureCircle USA Inc., Oak Brook, IL
- -Miller Coors / Crispin Cider, Colfax, CA
- -Denali Ingredients, New Berlin, WI
- -Griffith Labs, Chicago, IL
- -Rich Products, Niles, IL -Stratas Foods, Dallas, TX
- -Clos du Bois, Geyserville, CA
- -Cargill Foods, US Nationwide
- -Quaker Oats Company

Educational

- -Highland Community College, Freeport, IL
- -Northwest College, Powell, WY
- -Kansas State University, Manhattan, KS
- -University of Missouri, Columbia, MO
- -Southern Illinois University, Carbondale, IL
- -Wayne County Community College District, Belleville, MI
- -Claremont Colleges, Claremont, CA
- -Andrews University, Berrien Springs, MI
- -University of St. Francis, Joliet, IL
- -Leman Manhattan Prep School, New York, NY
- -Marquette University, Milwaukee, WI

Medical / Health Care

- -Los Alamos Medical Center, Los Alamos, NM
- -Stryker, Fremont, CA
- -Triclinic Labs, Lafayette, IN
- -North Western Memorial Hospital, Chicago, IL
- -North Star Medical Radioisotopes, Beloit, WI -Loreto Hospital, Chicago, IL
- -Ayass Lung Clinic & Sleep Center, Frisco, TX
- -Xenex Disinfection Services, San Antonio, TX
- -Fresenius Kabi, Lake Zurich, IL
- -Sage Products Inc., Cary, IL
- -Wheaton Franciscan Medical Center, Franklin, WI
- -Baxter Health Care, Los Angeles, CA
- -Care Fusion, Vernon Hills, IL
- -Mount Sinai School of Medicine, New York, NY

Industrial / Manufacturing

- -Ingersoll Rand (Trane Company), La Crosse, WI
- -Aerojet/Rocketdyne, Jonesborough, TN
- -IP Corp, St. Paul, MN
- -Lanco Pennland, Hancock, MD
- -Michelin, US Nationwide
- -Advanced Atomization Technologies, Clyde, NY
- -Enviroplas, Evansville, IN
- -Xerox, Webster, NY
- -Schneller, Kent, OH
- -ARC Products, Inc., Dallas, TX
- -Guardian Industries, Carleton, MI
- -KIK Products, Danville, IL
- -Bausch & Lomb, Clearwater, FL
- -Omni Logistics, Coppell, TX
- -Parkersburg Tile Co., Parkersburg, WV
- -Oberlin Filters, Waukesha, WI
- -Saint-Gobain, San Diego, CA
- -CertainTeed, Norwood, MA
- -The Honest Company, Playa Vista, CA
- -Whitehouse Laboratories, Whitehouse, NJ
- -Domtar Manufacturing, Bennettsville, SC
- -Federal-Mogul, Waupun, WI
- -Buckeye Florida Paper, Perry, FL
- -Miromatrix, Glencoe, MN

Chemical & Plastics

- -Oil-Dri America, Ochlocknee, GA
- -Chemtool, Division of Lubrizol, Rockton, IL
- -Mississippi Silicone, Burnsville, MS
- -Valspar, Wheeling, IL
- -Vantage Specialty Chemical, Gurnee, IL
- -Akzo Nobel, Reading, PA
- -AOC Resins, Collierville, TN
- -Enviroplas, Evansville, IN
- -Oakland Chemical, Rochester, MN
- -Petroleum Accessories, Houston, TX
- -Brenntag Great Lakes, LLC, Menomonee Falls, WI
- -Sigma Aldrich, Sheboygan Falls, WI
- -3M, Cordova, IL
- -Rust-Oleum Corp. Pleasant Prairie, WI

OUR CUSTOMERS -NATION AND WORLD WIDE!

Aggregate & Mining

- -Krete Industries, Brookfield, WI
- -Titan America, Medley, FL
- -Morton Salt, Elgin, IL & Port Canaveral, FL
- -Mohawk Asphalt, Skotia, NY
- -Boral, San Antonio, TX

Metallurgical

- -Fisher Barton, Watertown, WI
- -Lap Master International, LLC, Mount Prospect, IL
- -Woodward, Rockford, IL
- -Ingersoll Rand (Trane Company), La Crosse, WI
- -Hainbuch America, Germantown, WI
- -Heraeus Precious Metals, W. Conshohocken, PA

Water & Waste Water

- -Rock-Tenn, Delaware Water Gap, PA
- -Veolia Water, Yukon, OK
- Urbana, IL
 -North Liberty WTP, North Liberty, IN
- -Colorado Springs WTP, Colorado Springs, CO

-Fisher Scientific / IL State Water Survey,

- -Muscatine Water Pollution Control Plant, Muscatine, IA
- -Village of Tonica, Tonica, IL

International

Saudi Arabia

- -Gulf Keystone Petroleum, Kurdistan, Iraq -Hyundai Engineering & Construction,
- -Tatweer Mud Lab, Manama City, Bahrain
- -Diageo, Panama City, Panama
- -Diageo, St Croix, British V.I.
- -Western Arctic Research Center, Inuvik, Canada
- -Exxon Mobile, Papua New Guinea
- -Entomology Lab, Iquitos, Peru